


Silimed Tissue Expanders - Breast




Silimed Round

SMOOTH

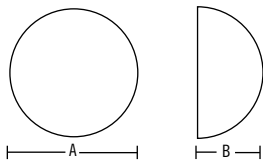

| REF | VOLUME ml | DIMENSIONS cm | | VALVE OPTIONS | | | | |
|-----|--------------|------------------|-----|---------------|---|---|--|--|
| | | A | B | | | | | |
| 800 | 50 | 6.2 | 3.6 | • | • | • | | |
| | 100 | 7.9 | 2.8 | • | • | • | | |
| | 200 | 9.6 | 4.1 | • | • | • | | |
| | 300 | 10.6 | 5.0 | • | • | • | | |
| | 400 | 11.2 | 5.5 | • | • | • | | |
| | 500 | 12.2 | 6.2 | • | • | • | | |
| | 600 | 13.2 | 6.5 | • | • | • | | |
| | 700 | 15.0 | 6.1 | • | • | • | | |
| | 800 | 15.2 | 6.4 | • | • | • | | |
| | 900 | 15.2 | 7.6 | • | • | • | | |
| | 1000 | 15.3 | 8.4 | • | • | • | | |

TEXTURED

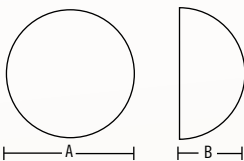


| REF | VOLUME ml | DIMENSIONS cm | | VALVE OPTIONS | | | | |
|-------|--------------|------------------|-----|---------------|--|--|--|---|
| | | A | B | | | | | |
| 20809 | 200 | 9.6 | 4.1 | | | | | • |
| | 300 | 10.6 | 5.0 | | | | | • |
| | 400 | 11.2 | 5.5 | | | | | • |
| | 500 | 12.2 | 6.2 | | | | | • |
| | 600 | 13.2 | 6.5 | | | | | • |
| | 700 | 15.0 | 6.1 | | | | | • |
| | 800 | 15.2 | 6.4 | | | | | • |
| | 900 | 15.2 | 7.6 | | | | | • |
| | 1000 | 15.3 | 8.4 | | | | | • |

Tissue Expander Valves


We recommend the use of a 21 gauge needle, or smaller.

VALVE OPTIONS:

- = Round (adult)
- = Round (child)
- = Luer Lock
- = Cone Shaped
- = Incorporated Magnetic

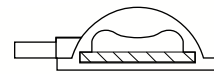

ROUND ADULT VALVE:
Base Diameter 34mm

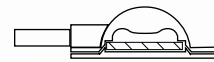

ROUND CHILD VALVE:
Base Diameter 22mm

Round Valves have a resistant radiopaque stainless steel disk to prevent the needle from penetrating the valve base. The needle must be introduced all the way to the bottom perpendicular to the valve.




LUER LOCK:

This luer accompanies intra-operative expanders and some expanders for slow expansion as an external valve.




CONE VALVE:

Dome Diameter 11mm and Base Diameter 22mm. Each contains a ceramic cup to prevent the needle from penetrating the valve base.

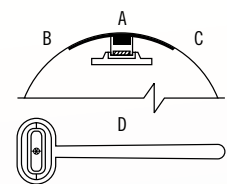

INCORPORATED MAGNETIC VALVE FOR SLOW EXPANSION WITH A MAGNETIC FINDER:


This valve has a Silicone Shell Septum (A), reinforcement to prevent folds (B), magnetic disk with stabilizer (C), and comes with a magnetic finder (D).

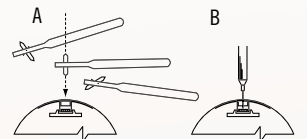

LOCATING THE VALVE:

(A) Marking the exact place for the injection. (B) In introducing the needle, make sure it penetrated the septum completely

IMPORTANT: Whenever necessary, we can supply any kind of expander with a Child or Adult Round Valve, Cone Valve or Luer Lock.

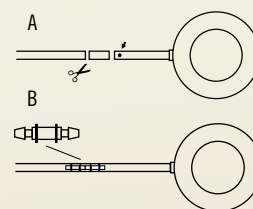

ASSEMBLY SEQUENCE:


(A) Cut the tube to the required length without going past the mark. (B) Place expander tube into the free end of the connective pin until it meets the other tube.



Federal (USA) law restricts this device to sale by or on the order of a physician.


Silimed Crescent

SMOOTH

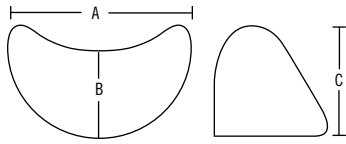

| REF | VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|-----|-----------|---------------|-----|-----|---------------|---|---|---|---|
| | | A | B | C | | | | | |
| 820 | 12 | 5.1 | 2.1 | 2.1 | • | • | • | | |
| | 25 | 6.0 | 2.7 | 2.6 | • | • | • | | |
| | 50 | 6.8 | 3.5 | 3.3 | • | • | • | | |
| | 100 | 8.1 | 4.3 | 3.5 | • | • | • | | |
| | 150 | 8.7 | 4.6 | 4.1 | • | • | • | | |
| | 200 | 9.3 | 4.9 | 4.7 | • | • | • | | |
| | 300 | 10.1 | 5.6 | 5.8 | • | • | • | | |
| | 400 | 11.4 | 6.2 | 6.5 | • | • | • | | |
| 700 | 13.4 | 7.8 | 8.8 | • | • | • | | | |
| 827 | 12 | 5.1 | 2.1 | 2.1 | | | | | • |
| | 25 | 6.0 | 2.7 | 2.6 | | | | | • |
| | 50 | 6.8 | 3.5 | 3.3 | | | | | • |
| | 100 | 8.1 | 4.3 | 3.5 | | | | | • |
| | 150 | 8.7 | 4.6 | 4.1 | | | | | • |
| | 200 | 9.3 | 4.9 | 4.7 | | | | | • |
| | 300 | 10.1 | 5.6 | 5.8 | | | | | • |
| | 400 | 11.4 | 6.2 | 6.5 | | | | | • |
| 700 | 13.4 | 7.8 | 8.8 | | | | | • | |
| 829 | 150 | 8.7 | 4.6 | 4.1 | | | | | • |
| | 200 | 9.3 | 4.9 | 4.7 | | | | | • |
| | 300 | 10.1 | 5.6 | 5.8 | | | | | • |
| | 400 | 11.4 | 6.2 | 6.5 | | | | | • |
| | 700 | 13.4 | 7.8 | 8.8 | | | | | • |

TEXTURED

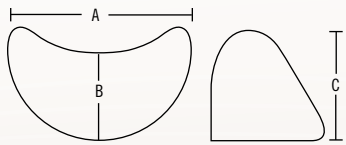

| REF | VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|-------|-----------|---------------|-----|-----|---------------|---|---|---|---|
| | | A | B | C | | | | | |
| 20820 | 12 | 5.1 | 2.1 | 2.1 | • | • | • | | |
| | 25 | 6.0 | 2.7 | 2.6 | • | • | • | | |
| | 50 | 6.8 | 3.5 | 3.3 | • | • | • | | |
| | 100 | 8.1 | 4.3 | 3.5 | • | • | • | | |
| | 150 | 8.7 | 4.6 | 4.1 | • | • | • | | |
| | 200 | 9.3 | 4.9 | 4.7 | • | • | • | | |
| | 300 | 10.1 | 5.6 | 5.8 | • | • | • | | |
| | 400 | 11.4 | 6.2 | 6.5 | • | • | • | | |
| 700 | 13.4 | 7.8 | 8.8 | • | • | • | | | |
| 20827 | 12 | 5.1 | 2.1 | 2.1 | | | | | • |
| | 25 | 6.0 | 2.7 | 2.6 | | | | | • |
| | 50 | 6.8 | 3.5 | 3.3 | | | | | • |
| | 100 | 8.1 | 4.3 | 3.5 | | | | | • |
| | 150 | 8.7 | 4.6 | 4.1 | | | | | • |
| | 200 | 9.3 | 4.9 | 4.7 | | | | | • |
| | 300 | 10.1 | 5.6 | 5.8 | | | | | • |
| | 400 | 11.4 | 6.2 | 6.5 | | | | | • |
| 700 | 13.4 | 7.8 | 8.8 | | | | | • | |
| 20829 | 150 | 8.7 | 4.6 | 4.1 | | | | | • |
| | 200 | 9.3 | 4.9 | 4.7 | | | | | • |
| | 300 | 10.1 | 5.6 | 5.8 | | | | | • |
| | 400 | 11.4 | 6.2 | 6.5 | | | | | • |
| | 700 | 13.4 | 7.8 | 8.8 | | | | | • |



Federal (USA) law restricts this device to sale by or on the order of a physician.




Silimed Tissue Expanders - Non-Breast



Silimed Auricular

SMOOTH

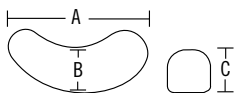


| REF | VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|-----|--------------|------------------|-----|-----|---------------|--|--|--|---|
| | | A | B | C | | | | | |
| 841 | 16 | 6.2 | 1.6 | 1.1 | | | | | • |
| | 40 D/E | 6.5 | 2.2 | 2.6 | | | | | • |
| | 55 D/E | 7.0 | 2.5 | 3.2 | | | | | • |

D=right, E=left

Silimed Calf


SMOOTH

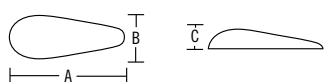

| REF | VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|-----|--------------|------------------|-----|-----|---------------|--|--|--|---|
| | | A | B | C | | | | | |
| 880 | 145 | 19.2 | 5.3 | 2.6 | | | | | • |
| | 185 | 22.9 | 5.7 | 2.6 | | | | | • |

Silimed Finger-Syndactylism

SMOOTH

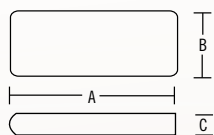

| REF | VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|-----|--------------|------------------|-----|-----|---------------|--|--|--|---|
| | | A | B | C | | | | | |
| 851 | 1 | 1.0 | 1.0 | 0.3 | | | | | • |
| | 2 | 2.0 | 1.0 | 0.3 | | | | | • |
| | 3 | 3.0 | 1.0 | 0.3 | | | | | • |
| | 4 | 4.0 | 1.0 | 0.3 | | | | | • |
| | 5 | 5.0 | 1.0 | 0.3 | | | | | • |



Federal (USA) law restricts this device to sale by or on the order of a physician.


*WR – When the desired expander does not require reinforcement, add the letters WR (without reinforcement) to the reference number. Ex: Rectangular expander, 160 ml, round adult valve, without reinforcement = 810-160WR

Tissue Expander Valves


We recommend the use of a 21 gauge needle, or smaller.

VALVE OPTIONS:

- = Round (adult)
- = Round (child)
- = Luer Lock
- = Cone Shaped
- = Incorporated Magnetic

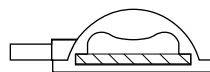

ROUND ADULT VALVE:
Base Diameter 34mm

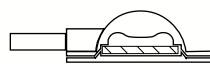


ROUND CHILD VALVE:
Base Diameter 22mm

Round Valves have a resistant radiopaque stainless steel disk to prevent the needle from penetrating the valve base. The needle must be introduced all the way to the bottom perpendicular to the valve.




LUER LOCK:

This luer accompanies intra-operative expanders and some expanders for slow expansion as an external valve.

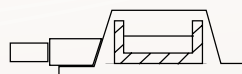


CONE VALVE:

Dome Diameter 11mm and Base Diameter 22mm. Each contains a ceramic cup to prevent the needle from penetrating the valve base.

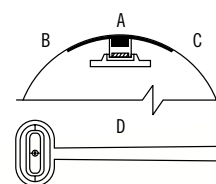

INCORPORATED MAGNETIC VALVE FOR SLOW EXPANSION WITH A MAGNETIC FINDER:


This valve has a Silicone Shell Septum (A), reinforcement to prevent folds (B), magnetic disk with stabilizer (C), and comes with a magnetic finder (D).

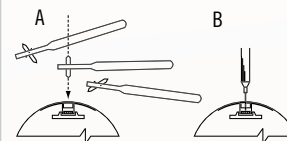

LOCATING THE VALVE:

(A) Marking the exact place for the injection. (B) In introducing the needle, make sure it penetrated the septum completely

IMPORTANT: Whenever necessary, we can supply any kind of expander with a Child or Adult Round Valve, Cone Valve or Luer Lock.

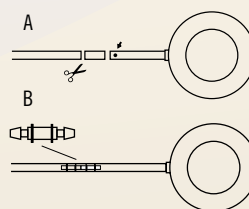


ASSEMBLY SEQUENCE:

(A) Cut the tube to the required length without going past the mark. (B) Place expander tube into the free end of the connective pin until it meets the other tube.

Silimed Longitudinally Curved


SMOOTH

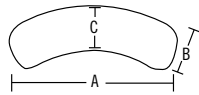

| REF | VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|-----|-----------|---------------|-----|-----|---------------|---|---|--|--|
| | | A | B | C | | | | | |
| 830 | 220 | 10.0 | 4.0 | 4.5 | • | • | • | | |
| | 480 | 20.7 | 5.0 | 4.8 | • | • | • | | |
| | 520 | 16.0 | 5.0 | 6.0 | • | • | • | | |

Silimed Oval/Round


TEXTURED

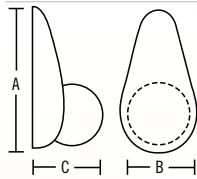

| REF | ANTERIOR / POSTERIOR VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|-----------|--------------------------------|---------------|------|-----|---------------|--|--|--|---|
| | | A | B | C | | | | | |
| 20739-320 | 140/180 | 15.0 | 10.5 | 5.6 | | | | | • |
| 20739-470 | 220/250 | 15.5 | 11.6 | 5.8 | | | | | • |
| 20739-600 | 300/300 | 16.3 | 12.3 | 6.5 | | | | | • |
| 20739-850 | 400/450 | 17.6 | 14.1 | 6.7 | | | | | • |


Silimed Rectangular

SMOOTH

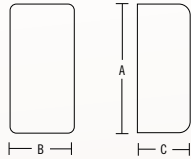

| REF | VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|------|-----------|---------------|-----|-----|---------------|---|---|--|---|
| | | A | B | C | | | | | |
| 810 | 30 | 5.1 | 2.5 | 2.5 | • | • | • | | |
| | 50 | 3.8 | 3.8 | 3.8 | • | • | • | | |
| | 100 | 6.0 | 4.6 | 5.3 | • | • | • | | |
| | 105 | 8.7 | 4.1 | 3.3 | • | • | • | | |
| | 160 | 11.0 | 5.1 | 3.8 | • | • | • | | |
| | 200 | 6.0 | 5.0 | 7.0 | • | • | • | | |
| | 205 | 9.5 | 5.5 | 4.0 | • | • | • | | |
| | 250 | 9.6 | 5.9 | 5.0 | • | • | • | | |
| | 320 | 14.1 | 6.2 | 4.7 | • | • | • | | |
| | 400 | 13.6 | 5.5 | 5.6 | • | • | • | | |
| | 500 | 15.8 | 6.9 | 5.1 | • | • | • | | |
| | 600 | 18.0 | 5.0 | 7.0 | • | • | • | | |
| | 640 | 18.0 | 8.0 | 5.0 | • | • | • | | |
| | 700 | 28.0 | 5.0 | 5.0 | • | • | • | | |
| | 720 | 15.2 | 7.3 | 6.3 | • | • | • | | |
| 800 | 18.5 | 8.7 | 5.7 | • | • | • | | | |
| 1000 | 20.2 | 8.0 | 7.0 | • | • | • | | | |
| 819 | 160 | 11.0 | 5.1 | 3.8 | | | | | • |
| | 200 | 6.0 | 5.0 | 7.0 | | | | | • |
| | 205 | 9.5 | 5.5 | 4.0 | | | | | • |
| | 250 | 9.6 | 5.9 | 5.0 | | | | | • |
| | 320 | 14.1 | 6.2 | 4.7 | | | | | • |
| | 400 | 13.6 | 5.5 | 5.6 | | | | | • |
| | 500 | 15.8 | 6.9 | 5.1 | | | | | • |
| | 600 | 18.0 | 5.0 | 7.0 | | | | | • |
| | 640 | 18.0 | 8.0 | 5.0 | | | | | • |
| | 700 | 28.0 | 5.0 | 5.0 | | | | | • |
| | 720 | 15.2 | 7.3 | 6.3 | | | | | • |
| | 800 | 18.5 | 8.7 | 5.7 | | | | | • |
| | 1000 | 20.2 | 8.0 | 7.0 | | | | | • |

TEXTURED

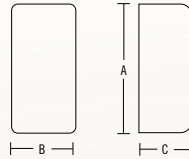

| REF | VOLUME ml | DIMENSIONS cm | | | VALVE OPTIONS | | | | |
|-------|-----------|---------------|-----|-----|---------------|---|---|--|---|
| | | A | B | C | | | | | |
| 20810 | 30 | 5.1 | 2.5 | 2.5 | • | • | • | | |
| | 50 | 3.8 | 3.8 | 3.8 | • | • | • | | |
| | 100 | 6.0 | 4.6 | 5.3 | • | • | • | | |
| | 105 | 8.7 | 4.1 | 3.3 | • | • | • | | |
| | 160 | 11.0 | 5.1 | 3.8 | • | • | • | | |
| | 200 | 6.0 | 5.0 | 7.0 | • | • | • | | |
| | 205 | 9.5 | 5.5 | 4.0 | • | • | • | | |
| | 250 | 9.6 | 5.9 | 5.0 | • | • | • | | |
| | 320 | 14.1 | 6.2 | 4.7 | • | • | • | | |
| | 400 | 13.6 | 5.5 | 5.6 | • | • | • | | |
| | 500 | 15.8 | 6.9 | 5.1 | • | • | • | | |
| | 600 | 18.0 | 5.0 | 7.0 | • | • | • | | |
| | 640 | 18.0 | 8.0 | 5.0 | • | • | • | | |
| | 700 | 28.0 | 5.0 | 5.0 | • | • | • | | |
| | 720 | 15.2 | 7.3 | 6.3 | • | • | • | | |
| 800 | 18.5 | 8.7 | 5.7 | • | • | • | | | |
| 1000 | 20.2 | 8.0 | 7.0 | • | • | • | | | |
| 20819 | 160 | 11.0 | 5.1 | 3.8 | | | | | • |
| | 200 | 6.0 | 5.0 | 7.0 | | | | | • |
| | 205 | 9.5 | 5.5 | 4.0 | | | | | • |
| | 250 | 9.6 | 5.9 | 5.0 | | | | | • |
| | 320 | 14.1 | 6.2 | 4.7 | | | | | • |
| | 400 | 13.6 | 5.5 | 5.6 | | | | | • |
| | 500 | 15.8 | 6.9 | 5.1 | | | | | • |
| | 600 | 18.0 | 5.0 | 7.0 | | | | | • |
| | 640 | 18.0 | 8.0 | 5.0 | | | | | • |
| | 700 | 28.0 | 5.0 | 5.0 | | | | | • |
| | 720 | 15.2 | 7.3 | 6.3 | | | | | • |
| | 800 | 18.5 | 8.7 | 5.7 | | | | | • |
| | 1000 | 20.2 | 8.0 | 7.0 | | | | | • |



PO Box 1490 | Santa Barbara, CA 93116-1490 | Tel: 805.562.3500 | Toll Free: 888.708.0808

Fax: 805.562.8401 | www.sientra.com

MDOC-DSNB Rev 001 9/2008